



3603 Ultra Slide

Surgical Table



Delivering 26" of Top Slide Without Any Compromises



Pairing the industry's longest top slide with substantial weight capacities, imaging access is only part of the 3603 Ultra Slide Series advantage. Combine the table's top slide capabilities with an 800 pound (363kg) articulation capacity, and the result is the best of both worlds: imaging and power. The table's impressive 26" (66 cm) of top slide accommodates the industry's most generous articulation capacity and provides superior imaging for patients of diverse sizes.

Imaging Efficiency


The 3603 Ultra Slide Series makes full-body imaging easier. At the touch of a button, patients can be effortlessly moved through the imaging window. The low-profile base stays clear of the imaging equipment, providing hassle-free C-arm access.

Surgical Positioning Flexibility

The removable back and leg section allow for easy and efficient accessory attachment. Enhance imaging capabilities with the optional 40" (102 cm) carbon fiber extension or attach the beach chair positioner for shoulder surgeries. With the back section removed, the beach chair can be fully integrated and powered by the table, providing one-touch setup.

Improved Patient Safety

Capable of 26" (66 cm) top slide at a full articulation capacity of 800 pounds (363 kg), the 3603 Ultra Slide Series will remain secure for patients of larger sizes. Its extensive capabilities not only make patient imaging simple, but also eliminate the need for risky repositioning.

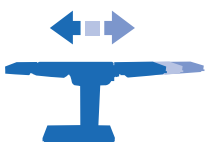

Market
Leading
Top Slide



Powerful
Articulation
Capacity



One-Touch
Imaging
Access



3603 Ultra Slide Series

26" (66 cm) Top Slide, General Purpose Table

Offering the industry's maximum top slide, the 3603 Ultra Slide Series takes the work out of imaging. At the click of a button, the patient can be moved throughout the expansive imaging window.



Differentiating Features

- 26" (66 cm) top slide
- 10" (25 cm) slide to head and 16 ½" (42 cm) slide to foot
- 1,000 lb (454kg) lift and 800 lb (363kg) articulation weight capacity
- Removable back and leg sections for quick accessory attachment and positioning
- Power beach chair positioner for one-touch set-up
- Durable stainless steel base is easy to clean
- Adjustable table height
 - 3603 20 ½" to 39 ½" (52 cm to 100 cm)
 - 3603H 23" to 41 ½" (58 cm to 105 cm)
- 35° lateral tilt
- 40° Trendelenburg and 30° reverse Trendelenburg




More Top Slide and Power in One Table

Able to support a 1,000 pound (454kg) lift and 800 pound (363kg) articulation weight capacity, the 3603 Ultra Slide Series provides 26" (66 cm) of top slide with greater power. Its superior design supports the full articulation capacity even during 26" (66 cm) of top slide – something the competition cannot do. As the table top slides toward the foot end, the head section can be conveniently detached and applied as a foot extension. With removable back and leg sections, the table offers quick attachment of accessories like the 40" (102 cm) carbon fiber extension and power beach chair positioner. Pre-programmed positioning options make setup even easier, allowing for one-touch maneuvers at any point throughout the procedure.

The 3603 Ultra Slide Series also brings confidence to the operating room. Recognized for their reliable performance, Skytron surgical tables are consistently rated as having a lower cost of ownership.

Lowest height on the market. There is the option of the 3603H, which offers a greater maximum height.

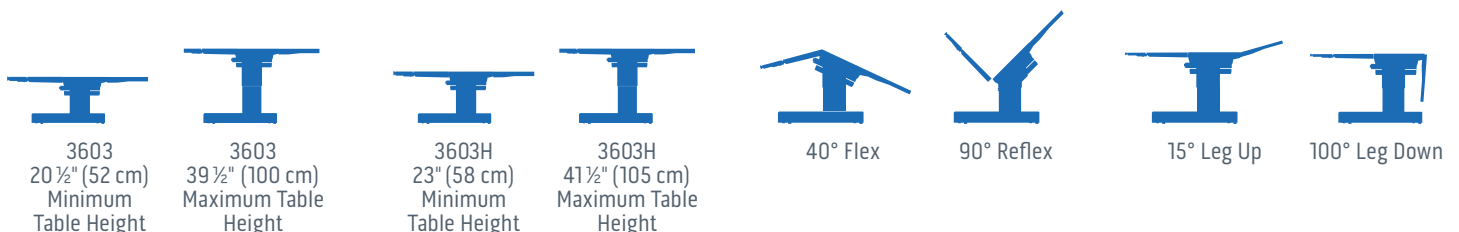



PSR - Patient Safe Return
Simultaneous return to level



Back lit back-up controls for emergency situations



New optional tunnel style
x-ray tops



Self-leveling brakes for stability
on uneven floors

One-Touch Pendant Control

- Simple one-touch operation
- Backlit control for high visibility in low lighting
- Audible alarm for kidney bridge, back up articulation and top slide obstruction
- With PSR (Patient Safe Return)
- Anatomically correct beach chair position
- Recessed buttons prevent accidental activation
- Durable and water resistant
- Can be used in conjunction with optional three function foot control



Optional Wireless Control



90° Back Up



40° Back Down



40° Trendelenburg



30° Reverse Trendelenburg



35° Tilt - Left



35° Tilt - Right

Positioning Guide




Upper Body Imaging



Lower Body Imaging



Upper Body Imaging with
40" (102 cm) Carbon Fiber Extension

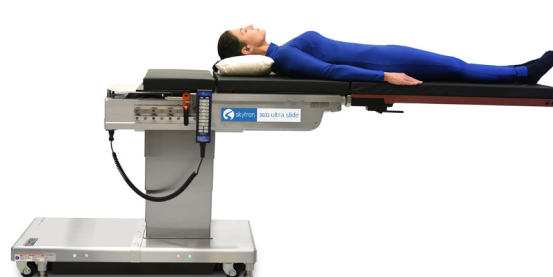

Lower Body Imaging with
40" (102 cm) Carbon Fiber Extension



Ophthalmic/ENT



Kidney/Thoracic

Bariatric Options



Bariatric



Abdominal



Neuro (Neck)



Neuro (Lumbar)



Shoulder Arthroscopy



Lap Nissen



Urology



Lumbar



Ophthalmic



Gastro/Intestinal

Skytron's Solutions



Sterile Processing

Processing Sinks and Stainless Offerings
Washers, Ultrasonics, and Cart Washers
Prep & Pack Tables & Sterilizers



Clinical

Surgical Tables
Booms & Lights
Full Offering of Stainless Products



Architectural & Environmental

Modular Walls
Configurable Window and Door Offerings
Built-in Cabinetry



Clinical Business Intelligence

Video Integration
Asset Tracking & Supply Chain Solutions
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service.

Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV3 03-24
PN L2-010-139