



Stellar XL

Surgical Light

 **SKYTRON**[®]
Reliability • Versatility • Cost Savings


High-Value Surgical Lighting with Exceptional Workflow Efficiency

Well suited for any surgical application, the Stellar XL light delivers adaptability to fit the needs of your facility. Skytron's unique reflector design technology provides a clear view with no obstructions during your procedure, allowing the surgeon to see precise detail.

Improved Visualization

Featuring a 27" range of total focusable depth, the Stellar XL's focal point can be easily fine-tuned with a twist of the center sterile handle – without moving the lighthead itself.

- Extend the working distance without moving the lighthead
- Reduce risk of contact with the lighthead or impeding line of sight to staff or displays
- Position the lighthead exactly where you need it with no drift

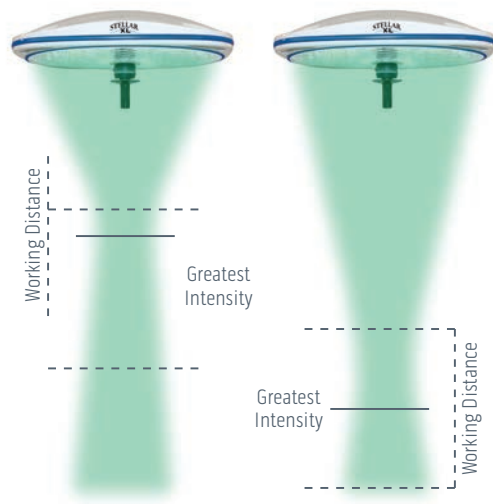

Eliminating Shadows During Procedures

Stellar XL provides the surgeon with a crisp, homogeneous spot and hidden LEDs that eliminate glare and shadows. Our improved thermal control allows for a cool, comfortable working environment. These features are accomplished with our unique Vertically Segmented Reflector Design Technology which:

- Utilizes high-intensity LEDs
- Redirects light rays into the surgical field using upward-facing LED technology
- Results in 2,000+ beams of light converging into a focusable column



Comfortable Work Environment

Stellar XL provides superior lighting and creates a comfortable working environment. The focusability of the light allows staff to position the lighthead farther away from the surgical site permitting more headroom for clinicians.



Eliminating Shadows



Effortless Operation



Comfortable Work Environment



Configurable For Your Unique Needs

Available in single, dual or triple lighthead configurations and fully compatible with Skytron booms, the Stellar XL is adaptable to many clinical settings.

Opportunities for Stellar XL

- Upgrading outdated technology
- Outpatient Surgery Centers
- General Operating Rooms
- C-Section Suites
- Labor and Delivery Suites
- Emergency Departments
- Intensive Care Units
- Procedure Rooms



Also available on a mobile stand

Stellar XL Details:

- Natural warm color temperature provides white light for high intensity applications
- 27" range of total focusable depth
- 24" diameter lighthead
- 360° infinite rotation on all axes permits unmatched flexibility and positioning
- Lightweight polymer housing permits easy maneuverability



Circle of Lighting Solutions



Aurora Four
True Focus Adjustment
Selectable Color Temperature



Lucina 4
LED Strobe-Guided Accuracy
Fully Recessed Design



Stellar XL
Crisp, Uniform Spotlighting
360° Rotation on All Axes



Spectra
Exceptional Exam Light
Excellent Shadow Reduction



Aurora Astro
Sterile Intensity Control
Premier Procedure Light



DoVera
Double Reflection Technology
Personalized Light Settings

Skytron is the Healthcare Efficiency Specialist, providing full-room solutions of capital equipment, architectural and real-time information systems for Medical, Surgical, Sterile Processing, and Infection Prevention. Our solutions enhance the utilization of people, facilities, and capital because they are designed with the user in mind and have a low, long-term cost of ownership.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV2 04-21
PN L2-010-124