



Integrity 215

Steam Sterilizers

Specifically engineered to offer excellent throughput in a space-saving footprint.

Integrity 215 exceeds all standards for capacity—efficient enough for immediate use, productive enough for high-demand environments, and compact enough for most replacement applications. Building on Skytron's history of quality and innovation, the 215 and 215SG are part of a family of sterilizers designed to meet the wide-ranging needs of today's healthcare environments. Whether you're outfitting a new facility, renovating an older one, or replacing outdated technology, Integrity sterilizers offer industry-changing efficiency, capacity, ease of use, and use of space.

Twice the capacity of similar sized sterilizers

Processes up to six 25-pound trays.

Easy to use, program, and check cycle status at a glance

Simple touchscreen controls and a large countdown timer and progress bar make daily operation smoother.

Save on energy, water, and chemicals

Integrated water recirculation reduces water usage by up to 67%*, without affecting cycle time or increasing footprint.

Long life and trouble-free operation

The Integrity Series offers fully jacketed stainless steel chambers, chain-driven electric-powered doors, and stainless steel and copper piping standard.

Superb diagnostics and support

System-wide parts monitoring, on-screen diagnostics, local technicians, and 24/7/365 call center ensure the reliable operation of your equipment.

Plug-and-play installation


Simplifies replacement of older technology.



Integrity 215


Steam Sterilizers

24 Preset and custom cycles

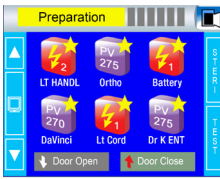

Processing cycles

- Immediate use 1 (3 minute)
- Immediate use 2 (4 minute)
- Express prevacuum (immediate use 3)
- Prevacuum 270
- Prevacuum 275
- Gravity 1
- Gravity 2
- Liquid cycle (optional)
- Preset cycles validated to AAMI ST-8:2008

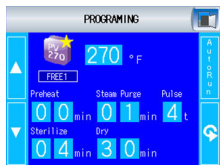


Test cycles

- Warm up
- BD
- Cycle
- Leak

Custom cycles

12 custom cycles modifiable to your unique needs

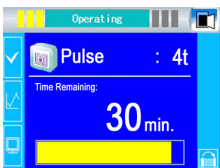

Sterilize quick turnaround devices or even complete orthopedic sets with simple touchscreen control

- Spacious 21½" x 21½" x 38" sterilization chamber
- Processes up to six 25-pound trays (150 pounds per cycle)
- Integrated water recirculation system
- 5.7" LCD full-color touchscreen panel
- Countdown timer and progression bar visible from more than 30'
- 12 standard and 12 fully customizable cycles
- Electric-powered sliding door
- Recessed foot-pedal and touchscreen door control
- Stainless steel and copper piping
- Dual-walled, fully jacketed stainless steel pressure vessel
- Programmed automatic cycle start
- Energy saving daily startup/shutdown timer
- Electronic recording of cycle report data
- Integrated ink-to-paper printer for hard copy cycle reports
- Password protected access to supervisor and service modes
- Compliant to AAMI-ST8:2008 and AAMI-ST79:2010 standards

Options

- Integrated steam generator (SG models)
- Load/unload cart and carriage
- Recessed mounting or freestanding cabinet.

Order Configurations

Integrity Series	Door Type	Heating	Control Power	Voltage (Opt)	Mounting	Equipment #
Integrity™ 215	Single	House Steam	120V, 60Hz, 10A	NA	Recessed	ASSV-AB09-120R
					Freestanding	ASSV-AB09-120F
Steam Generator		208V, 3pH, 60Hz, 75A		Recessed	ASSV-AB09-208R	
				Freestanding	ASSV-AB09-208F	
		480V, 3pH, 60Hz, 33A		Recessed	ASSV-AB09-480R	
				Freestanding	ASSV-AB09-480F	

*Water usage is based on comparisons of fully loaded sterilizers with similar-size chambers (20" x 20"), using a 270°F prevacuum cycle with 4-minute sterilize and 20- or 30-minute dry time. Testing showed that the Integrity 215 used substantially less water than the similarly sized STERIS® AMSCO® 400 (35% less) and STERIS CENTURY® (43% less).

Skytron is the Healthcare Efficiency Specialist, providing full-room solutions of capital equipment, architectural and real-time information systems for Medical, Surgical, Sterile Processing, and Infection Prevention. Our solutions enhance the utilization of people, facilities, and capital because they are designed with the user in mind and have a low, long-term cost of ownership.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
 Grand Rapids, MI 49512
 1.800.SKYTRON (759-8766)
 email: info@skytron.com
www.skytron.com



REV3 05-19
 PN L2-010-99