



3008 Clarity

Surgical Table

Simplify full body imaging and patient positioning



The 3008 Clarity brings Skytron quality and versatility to an imaging surgical table. With a radiolucent 83" carbon fiber table top and floating top functionality, the 3008 offers precise and infinite control of the top's movement. The Clarity's wide range of motion and radiolucency make it a versatile and adaptive table for a variety of specialties and case types.

Imaging Access

The 3008 Clarity is equipped with an 83" carbon fiber top capable of sliding in any direction. Fully radiolucent, the carbon fiber construction allows clinicians to use the lowest doses for imaging and still maintain excellent image clarity. The 3008 boasts an impressive 65" unobstructed imaging window.

Floating Top

The float function is engaged with a single button, allowing the user to reposition the top in any direction with ease even with a 500lb patient load.

Hand Control

A specially designed, comfort grip hand control is all that is needed to articulate the 3008 Clarity floating top. The float function is engaged with a single button, allowing the user to push or pull the top in any direction with ease.




Multiple
Direction
Top-Slide



Powerful
Articulation
Capacity



One-Button
Float Top
Control

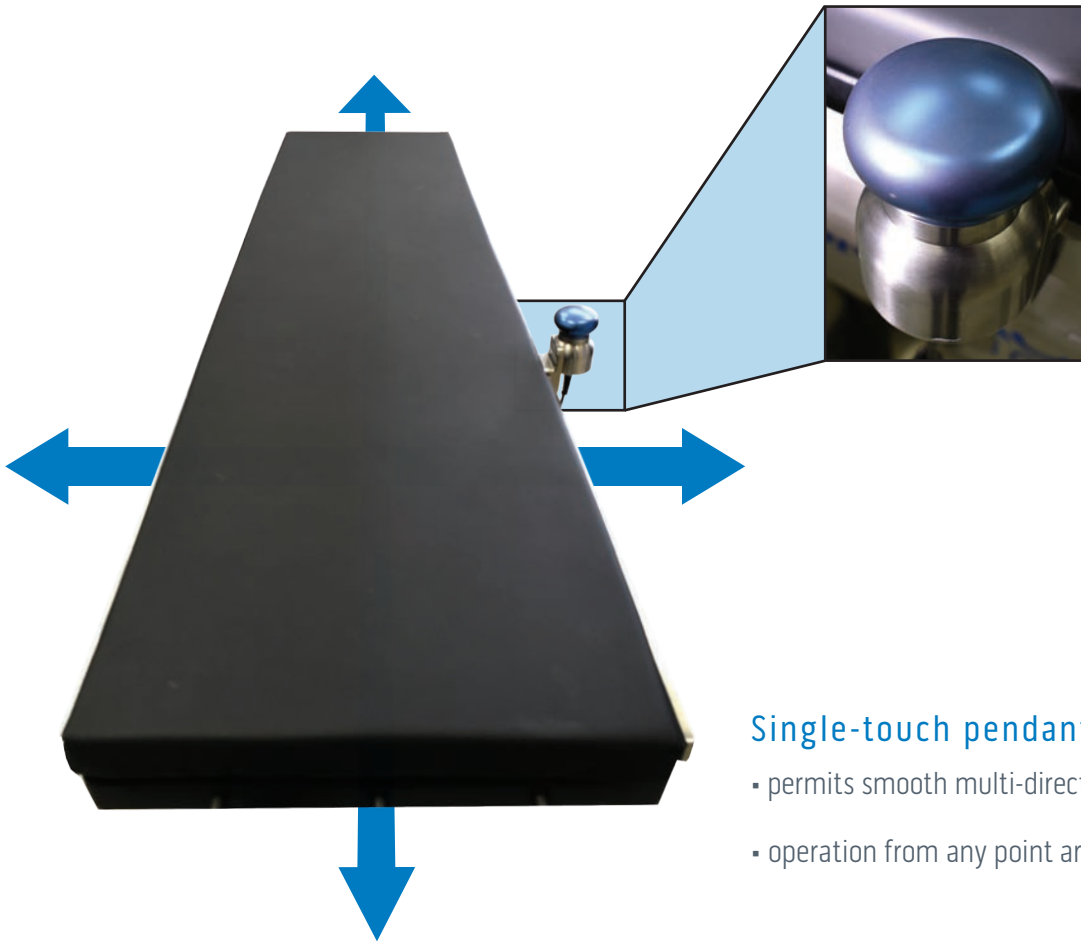

Single-touch pendant and free float controls

- permits smooth multi-directional top slide
- operation from any point around the table's perimeter

500 pound weight capacity includes lift and top slide

20° trend and tilt

Providing 20° of trend and tilt, the 3008 Clarity takes the work out of imaging setup. Radiolucent table top creates a minimum 65" of imaging length with added length possible with optional extensions. Tilt and trend support a 300 pound capacity through the entire range of motion.



Differentiating Features

- Radiolucent top for low dose images with excellent clarity
- Fluoroscopy made easy with lateral and longitudinal movement of floating top
- One button float top control
- 83" Carbon-fiber top for full body imaging
- 65" of unobstructed imaging length
- One touch return to level
- 500 lb lift and top slide capacity



Circle of Surgical Table Solutions



3603 Ultra Slide

High Performance
Removable Back
26½" Top Slide



6702 Hercules

Rotating Top &
Removable Back



3503 EZ Slide

High Performance
Top Slide



GS70 Salus

The Premier Solution
Smart, Safe, Simple



6302 Elite

Rotating Top



3008 Clarity

Premier Imaging Table



1602 Essentia

General Surgery Table

Skytron is the Healthcare Efficiency Specialist, providing full-room solutions of capital equipment, architectural and real-time information systems for Medical, Surgical, Sterile Processing, and Infection Prevention. Our solutions enhance the utilization of people, facilities, and capital because they are designed with the user in mind and have a low, long-term cost of ownership.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REVO 9-18
PN L2-010-156