



# DoVera DE12

Exam Light



## Enhance Your Performance with DoVera DE12


The DoVera DE12 is a cutting-edge exam light with fixed focus technology, designed to provide a shadow-free and glare-free environment. It features 12 elliptical parabolas that evenly distribute LED lights to eliminate any shadow formation. The light operates at a low temperature and provides color correct light at 4,350K, making it versatile for various procedures.

The DoVera DE12 is a highly advanced exam light that offers several key benefits:

- Slim profile and lightweight material for ease of use and storage
- Simple and intuitive intensity controls at your fingertips
- Easy to use side handles for quick positioning
- Sterilizable handle for positioning control during procedures
- Color correct light provides accurate visibility



Durable Construction

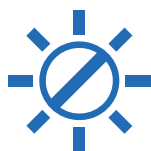

Reduced Eye Glare



Comfortable  
Side Handles

# Technical Specifications

Light intensity at 1 meter	119,330 lux
Lighthead diameter	16 inches
Color Temperature	4,350K
Color rendering index (CRI)	92 Ra
Spot Size	6.2 inches
Average LED life	60,000 hours
Control of illuminance	20% - 100%

## Several Mounting Options for Different Applications



### Dual Ceiling-Mounted

Dual ceiling mounted LED exam lights with 69" reach, 335° rotation of the upper radial arm and 360° rotation of the lower radial arm with 300° pitch and 180° roll rotation of each lighthead.



### Single Ceiling-Mounted

Single ceiling mounted LED exam lights with 66" reach, 360° arm rotation, 300° pitch, and 180° roll rotation of the lighthead.



### Wall Mounted

Wall mounted DE12 LED exam light with 66" reach, 180° rotation of the radial arm with 300° pitch and 180° roll rotation of the lighthead.



### Stand Light

DE12 LED exam light on mobile floor stand with 34" reach, 30° left and right rotation of the spring arm with 300° pitch and 180° roll rotation of the lighthead.

# Skytron's Solutions



## Sterile Processing

Processing Sinks and Stainless Offerings  
Washers, Ultrasonics, and Cart Washers  
Prep & Pack Tables & Sterilizers



## Clinical

Surgical Tables  
Booms & Lights  
Full Offering of Stainless Products



## Architectural & Environmental

Modular Walls  
Configurable Window and Door Offerings  
Built-in Cabinetry



## Clinical Business Intelligence

Video Integration  
Asset Tracking & Supply Chain Solutions  
Surgical Displays

Skytron's mission is to provide safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service. Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize unwavering transparency & accountability in our actions.

**To learn more, visit our website at [www.skytron.com](http://www.skytron.com)**



5085 Corporate Exchange Blvd. SE  
Grand Rapids, MI 49512  
1.800.SKYTRON (759-8766)  
email: [info@skytron.com](mailto:info@skytron.com)  
[www.skytron.com](http://www.skytron.com)



REV2 11-24  
PN L2-010-231