




DoVera
Surgical Light

 **SKYTRON**[®]
Reliability • Versatility • Cost Savings



Suited for any surgical application, Skytron's DoVera DM520 provides superior illumination while eliminating glare.

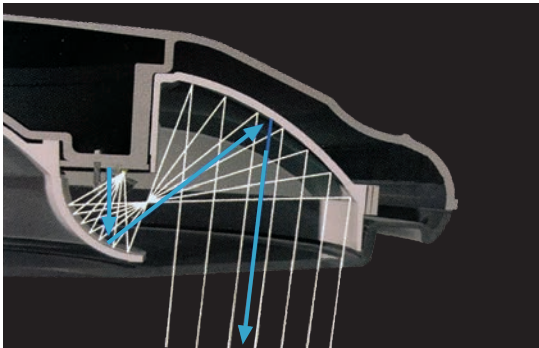


Double Reflection/No Glare

DoVera's unique Double Reflection technology eliminates glare at the lighthouse and provides uniform light across the surgical field. High intensity LED light beams project onto a segmented mirror, which reflects the light onto a smooth second mirror and subsequently to the surgical field. Double Reflection eliminates glare produced by traditional direct lighting technologies.




Personalized Lighting


Spot size and Color Temperature — Capacitive touch control at the yoke allow for easy, personalized set up of the light. Three spot sizes provide the right amount of light for every procedure while the seven color temperature options create a range of options from soft to bright white light. The surgical team can adjust the intensity to their preference throughout the procedure via the sterile center handle. Optional wall control available for adjustments outside of the field.

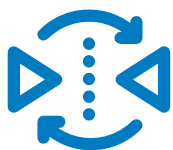

Double reflection technology

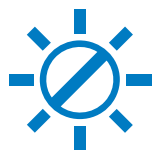

No glare

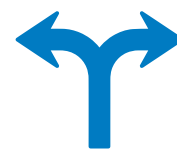

Personalized lighting



Intuitive controls on the yoke provide quick personalization of spot size, color temperature, and intensity.

Courtesy Light

Safety lighting for darkened surgical suites




Single Yoke

Easy positioning for low ceiling heights



Optional Wall Control

Non-sterile option for lighting adjustments

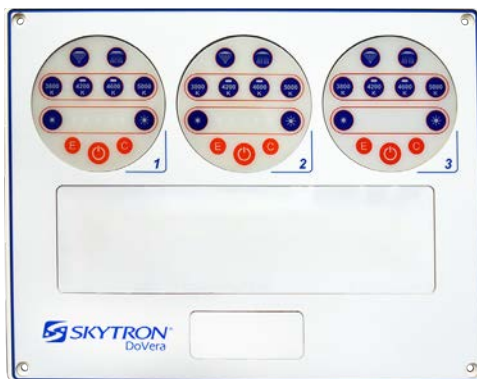

Double Yoke

Maximum positioning versatility





Circle of Lighting Solutions



Aurora Four
True Focus Adjustment
Selectable Color Temperature



Lucina 4
LED Strobe-Guided Accuracy
Fully Recessed Design



Spectra
Exceptional Exam Light
Excellent Shadow Reduction



DoVera
Double Reflection Technology
Personalized Light Settings



Stellar XL
Crisp, Uniform Spotlighting
360° Rotation on All Axes



Aurora Astro
Sterile Intensity Control
Premier Procedure Light

Skytron is the Healthcare Efficiency Specialist, providing full-room solutions of capital equipment, architectural and real-time information systems for Medical, Surgical, Sterile Processing, and Infection Prevention. Our solutions enhance the utilization of people, facilities, and capital because they are designed with the user in mind and have a low, long-term cost of ownership.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REVO 11-18
PN L2-010-153