



# Flex

Ultrasonic Washer-Disinfecter



## Clean Cannulated & Non-Cannulated in the Same Load

The Skytron Flex is an ultrasonic washer-disinfector designed specifically for orthopedic, robotic, cannulated, and non-cannulated surgical instruments.

### Cleaner


- Eliminates risks from manual cleaning
- Independent lab tested and validated
- Irrigates with fresh, clean water and detergent
- Supports compliance to AAMI ST79 Standards

### Faster

- Triple the capacity of standard ultrasonic washer-disinfectors
- Multi-step cycle cleans and disinfects in less than 40 minutes
- Non-critical water rinse cycles available to shorten cycle time
- Streamlined cleaning process for efficient results

### Smarter

- One Skytron Flex is more cost-efficient than multiple ultrasonics and washer-disinfectors
- Simplify operations with sonic only and sonic irrigation cycles
- Innovative technology drives down costs

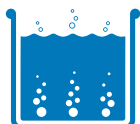


Ultrasonic  
Soak



Clean Water  
& Detergent  
Flush



Ozone  
Final  
Rinse

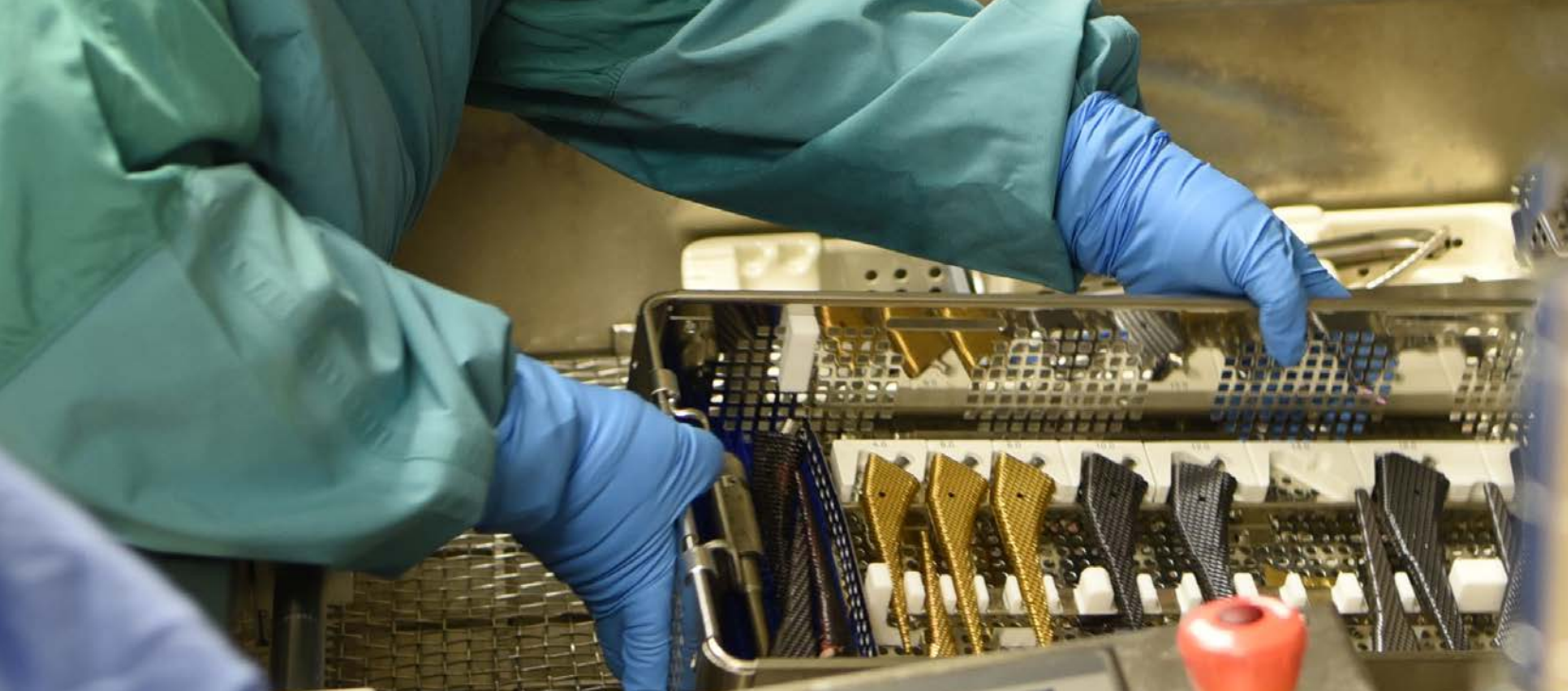

125 pound weight capacity includes up to six trays per cycle



### Differentiating Features

- Unsurpassed throughput
- Easy to use
- No proprietary trays-use what you currently have
- Fresh water used every cycle
- Cleans orthopedic, robotic, cannulated and non-cannulated instruments simultaneously

### Specifications

Dimensions	52" W x 41" D x 42" H
Weight Capacity	125 lbs
HF Output Power	2,000 watts
Frequency	40 kHz with a sweep of $\pm 2$ kHz
Washer Supply	480V/3 phase 50A, 60 Hz
Watts per Gallon	48 Watts/Gallon
Final Rinse	Critical Water Final Rinse



# Circle of Sterile Processing Solutions



## Integrity Sterilizers

Increased Efficiency  
Smaller Footprint



## Processing Sinks

Ergonomically Friendly Adjustable Height  
Customizable to Fit Specific Needs



## Flex Ultrasonics

Clean All Instruments in One Cycle  
Thermal Disinfection



## ErgoStat Prep & Pack Tables

Alleviate Physical Stress  
Height Adjustable Work Station



## Cart Washer-Disinfectors

Four Sizes to Choose From  
Quiet Operation



## Sterling Washer-Disinfectors

Water & Chemical Conservation  
Medical Device IFU Compliance

Skytron is the Healthcare Efficiency Specialist, providing full-room solutions of capital equipment, architectural and real-time information systems for Medical, Surgical, Sterile Processing, and Infection Prevention. Our solutions enhance the utilization of people, facilities, and capital because they are designed with the user in mind and have a low, long-term cost of ownership.

To learn more, visit our website at [www.skytron.com](http://www.skytron.com)



5085 Corporate Exchange Blvd. SE  
Grand Rapids, MI 49512  
1.800.SKYTRON (759-8766)  
email: [info@skytron.com](mailto:info@skytron.com)  
[www.skytron.com](http://www.skytron.com)



REV2 7-22  
PN L2-010-150