



Flex

Ultrasonic Washer-Disinfecter





Clean Cannulated & Non-Cannulated in the Same Load

The Skytron Flex is an ultrasonic washer-disinfector designed specifically for orthopedic, robotic, cannulated, and non-cannulated surgical instruments.

Reduce Labor and Error


- Eliminates staff risks from manual cleaning
- Independently lab tested and validated
- Irrigates with fresh, clean water and detergent
- Meets AAMI Standards
- Utilize your current baskets and vendor trays

Ease of Use

- Two thermal disinfection cycles - Cannulated and Non-Cannulated
- Reprocess mix loads in the same cycle
- Meets Intuitive Surgical Inc.'s recommended parameters for ultrasonic baths
- Built-in filtration system provides critical water final rinse

Efficient and Cost Effective

- Triple the capacity of standard ultrasonic washer-disinfectors
- Multi-step cycle cleans and disinfects in less than 35 minutes
- With its thermal disinfection rinse, Flex can be used as a backup washer-disinfector
- Reprocess 24 da Vinci instruments simultaneously

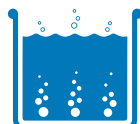

Ultrasonic
Soak



Clean Water
& Detergent
Flush



Ozone
Final
Rinse



125 Pound Weight Capacity Includes Up to Six Trays per Cycle



Differentiating Features

- Unsurpassed throughput
- Easy to operate
- No proprietary trays
- Fresh water used every cycle
- Cleans orthopedic, robotic, cannulated, and non-cannulated instruments simultaneously

Specifications

Dimensions	52" W x 41" D x 42" H
Basin Dimensions	40" W x 29" D x 10" H
Weight Capacity	125 lbs
HF Output Power	2,000 watts
Frequency	40 kHz with a sweep of ± 2 kHz
Washer Supply	480V/3 phase 50A, 60 Hz
Watts per Gallon	48 Watts/Gallon
Final Rinse	Critical Water Final Rinse

Skytron's Solutions



Sterile Processing

Processing Sinks and Stainless Offerings
Washers, Ultrasonics, and Cart Washers
Prep & Pack Tables & Sterilizers



Clinical

Surgical Tables
Booms & Lights
Full Offering of Stainless Products



Architectural & Environmental

Modular Walls
Custom Window and Door Offering
UVC Disinfection Products



Clinical Business Intelligence

Video Integration
Asset Tracking & Supply Chain Solutions
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service.

Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV4 1-23
PN L2-010-150