



Flex CT6

Countertop Ultrasonic Irrigator



The Solution Your SPD Deserves

With Sterile Processing Departments facing increasing pressure to reduce costs and improve efficiency, the CT6 assists in achieving these goals by providing a reliable, effective, and affordable way to clean surgical instruments.

Compact, Efficient, and Eco-Friendly


The CT6 is tailor-made for space-constrained facilities and consumes fewer utilities, contributing to both environmental sustainability and cost savings.

Patented Cleaning Excellence

By use of a patented suction and ultrasonic cleaning technique, the CT6 is meticulously designed to eliminate bioburden, debris, and contaminants from your surgical instruments. You can trust the CT6 to deliver impeccable cleanliness and safety.

Unparalleled Cleaning Standards

The CT6 exceeds the AAMI TIR 30 standard for cleaning effectiveness. When it comes to ensuring surgical instruments are at their absolute best quality, the CT6 delivers.

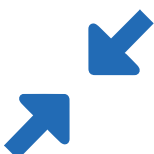

Small
Environmental
Footprint



Clean Water
& Detergent
Flush



Affordable
Ultrasonic
Solution



Where Innovation Meets Cleaning Excellence

Skytron's Flex CT6 Ultrasonic is more than just a device; it's a commitment to a safer, cleaner, and more efficient healthcare environment. It's time to revolutionize your sterile processing procedures and elevate your standards to new heights.

Similar to the Flex, the CT6 is engineered with versatility and throughput in mind. It handles various instrument types, from cannulated to non-cannulated, and even robotic instruments. Our patented dual hook-up method allows simultaneous cleaning of both ends, inside and out, ensuring a thorough, time-saving process. Plus, the CT6 can process up to 6 cannulated instruments at once, with each instrument having its own dedicated suction chamber, filter, and pump.



Specifications

| | |
|-------------------------|--|
| Dimensions | 33.5" W x 15.03" D x 16.63" H |
| Basin Dimensions | 30" W x 8" D x 6" H |
| Weight Capacity | 20 lbs |
| HF Output Power | 40 kHz |
| Frequency | 40 kHz |
| Washer Supply | 110/120 VAC 3 Cycle 50-60 Hz @ 9 Amps Max |
| Watts per Gallon | 80 Watts |

Skytron's Solutions



Sterile Processing

Processing Sinks and Stainless Offerings
Washers, Ultrasonics, and Cart Washers
Prep & Pack Tables & Sterilizers



Clinical

Surgical Tables
Booms & Lights
Full Offering of Stainless Products



Architectural & Environmental

Modular Walls
Window and Door Offerings
UVC Disinfection Products



Clinical Business Intelligence

Video Integration
Asset Tracking & Supply Chain Solutions
Surgical Displays

Skytron provides safe and efficient healthcare solutions delivered through a local representative network that values long-term relationships and comprehensive service.

Focused on providing the highest quality solutions while delivering the lowest cost of ownership, we prioritize transparency and accountability in our actions.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV1 11-23
PN L2-010-XXX