



Stellar XL

Surgical Light

 **SKYTRON**[®]
Reliability • Versatility • Cost Savings


High value surgical lighting with high-tech capabilities.

Well-suited for any surgical application, the Stellar XL light delivers advanced focusability and optimal shadow control to the operating room. The Stellar XL offers the ideal balance of features and value – whether updating outdated technology or outfitting a new space.

Handle Controlled Focus Capability

Featuring a 27" range of total focusable depth, the Stellar XL's focal point can be easily fine-tuned with the center sterile handle – without moving the lighthead itself.

- Set and forget the lighthead at a specific height
- Extend the working distance without moving the lighthead
- Reduce risk of contact with the lighthead or impeding line of sight to staff or displays

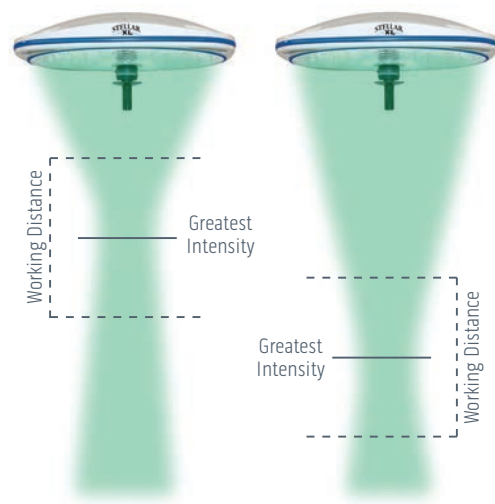


Vertically Segmented Reflector Design (VSRD)

Each of the high-intensity LED is suspended upward-facing in its cradle, allowing the vertically segmented reflector (VSRD) to redirect overlapping rays down to the surgical field. The result is 2,000+ beams of light that converge into a continuous column, leading to:

- A crisp, homogeneous spot
- Excellent shadow reduction
- Hidden LEDs that eliminate glare
- Improved thermal control
- A cool, comfortable working environment



True Focus Adjustment

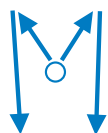

Vertically Segmented Reflector Design



Suspended LED System



Ideal For Any Application

Available in single, dual or triple lighthead configurations and fully compatible with Skytron booms, the Stellar XL is adaptable to many clinical settings.

Perfect for

- Upgrading outdated technology
- Outpatient Surgery Centers
- General Operating Rooms
- C-Section Suites
- Labor and Delivery Suites
- Emergency Departments
- Intensive Care Units
- Procedure Rooms




Also available on a mobile stand

Differentiating Features

- Natural warm color temperature provides white light for high intensity applications
- 27" range of total focusable depth
- 24" diameter lighthead
- 360° rotation on all axes permits unmatched flexibility and positioning
- Lightweight polymer housing permits easy maneuverability



Circle of Lighting Solutions

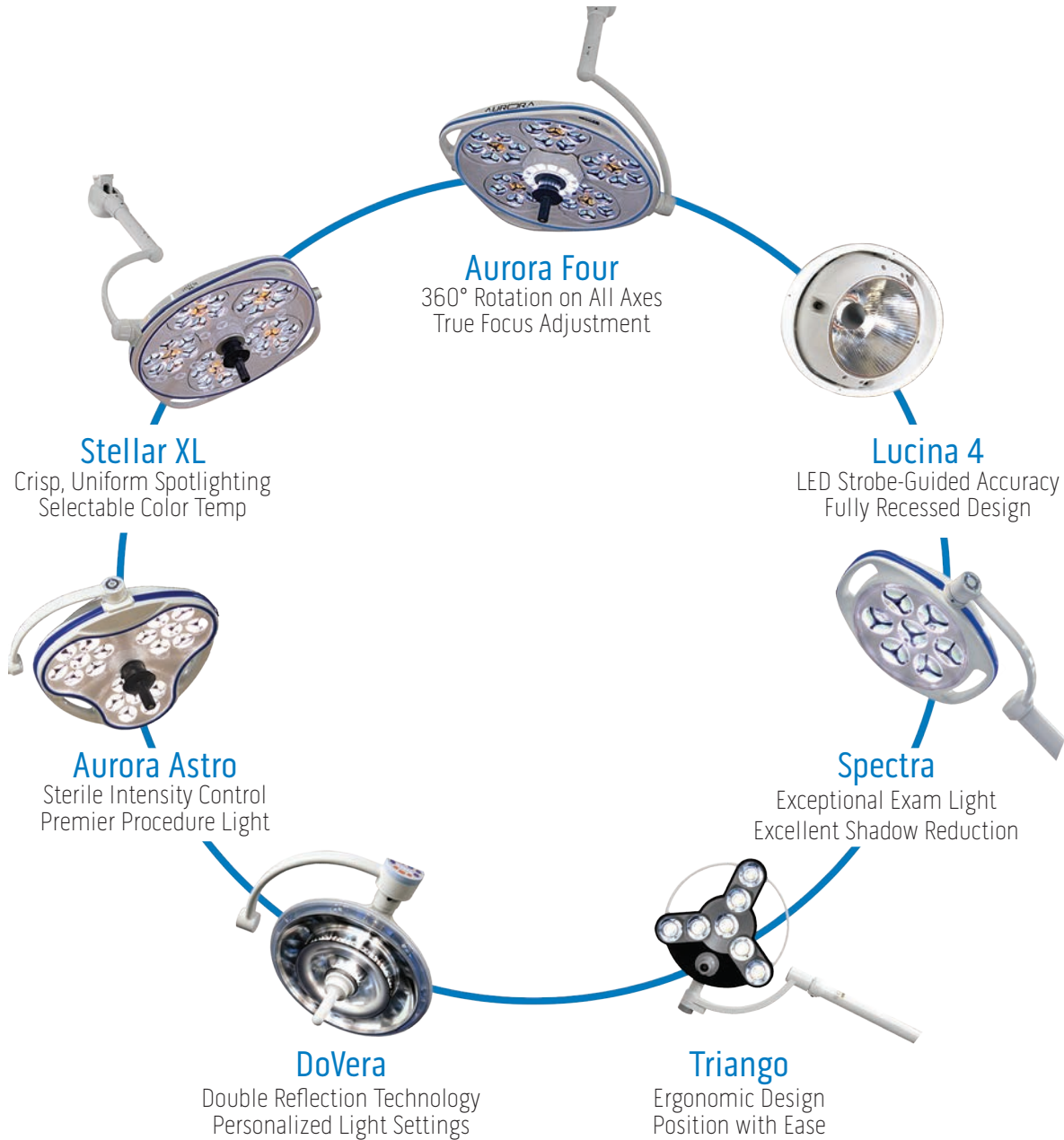

Skytron is the Healthcare Efficiency Specialist, providing full-room solutions of capital equipment, architectural and real-time information systems for Medical, Surgical, Sterile Processing, and Infection Prevention. Our solutions enhance the utilization of people, facilities, and capital because they are designed with the user in mind and have a low, long-term cost of ownership.

To learn more, visit our website at www.skytron.com



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV1 7-18
PN L2-010-124