



Sterling P1500

Washer-Disinfector Accessory Guide

Setting the standard for sterile processing...

The Sterling P1500 Washer-Disinfector streamlines and automates the decontamination process – bringing unsurpassed efficiency and flexibility to the Sterile Processing Department.

The P1500 offers:

Less Water Consumption - Averaging just 5.8 gallons (22L) per water fill, using up to 40% less water than its competition

Completely Open Cycle Architecture - 40 fully customized cycles, each with 10 available wash phases, providing the flexibility required to meet any thermal disinfection need

Process More in a Smaller Footprint- The industry's smallest overall design, with the deepest depth – provides the best throughput per cycle in the smallest footprint

Manual Transport Cart

Manual loading cart transports Sterling P1500 wash racks in all configurations. Easily maneuverable around all environments.

- Fits all Skytron Sterling P1500 wash rack configurations
- Built in drip tray and collector for excess water
- Easy to use coupling system for connection to the Sterling P1500 washer
- Cart dimensions are 25.8" W x 29.9"H x 32.3"L
(65.7 cm W x 76 cm H x 82.1 cm L)

Item #P1500-MT



Automatic Load/Unload

The Automatic loading system for the Sterling P1500 can be connected to the clean and unclean side of the machine.

- Fits all Skytron Sterling P1500 wash rack configurations
- Built in drip tray and collector for excess water
- Easy to use coupling system for connection to the Sterling P1500 washer
- Minimize heavy handling of instruments

Item #P1500-LOAD1

Item #P1500-UNLD1



Wash Racks

Flexible Rack

Two Levels

A flexible rack for multiple purposes. The flexible two-level rack can be combined with inserts for OR shoes, washing bowls and kidney bowls.

- Wash in up to 2 levels
- Includes grids and detachable washing arm
- Can be equipped with P1500-FLX-BWL, P1500-FLX-KB, P1500-B4L-BWL, and P1500-FLX-SH
- Maximum height lower level of 12.7 inches (32.4 cm) and 6.3 inches (15.9 cm) on the upper level

Item #P1500-FLX2



Flexible Rack

Four Levels

Universal 4 level rack suitable for various OR equipment from orthopedic surgeries - like hip, foot, knee, hand, shoulder or other surgeries that involve numerous instrument sets.

- Wash in up to 4 levels
- Includes grids and detachable washing arms
- Maximum height of 3.2 inches (8.2 cm) on the top level and 3.2 inches (8.2 cm) on each lower level

Item #P1500-FLX4



Flexible Rack

Five Levels

Universal 5 level rack suitable for various OR equipment from orthopedic surgeries - hip, foot, knee, hand, shoulder or other surgeries that involve many instruments and parts. Levels can be pulled out to make room for voluminous equipment, washing arms can be removed.

- Wash in up to 5 levels
- Includes grids and detachable washing arms
- Can wash up to 15 DIN instrument baskets per load
- Maximum height of 2.5 inches (6.4 cm) between levels

Item #P1500-FLX5



Wash Racks


Fixed Rack

Five Levels

Ideal solution when large capacity processing is a request. Uses five washing levels. Suitable for all kinds of OR-packages, like orthopedic surgery packages. Space for SPRI-baskets in various sizes.

- Washes up to 15 DIN instrument baskets per load.
- The loading space in each level is up to 2.5 inches (6.3 cm).

Item #P1500-FIX5

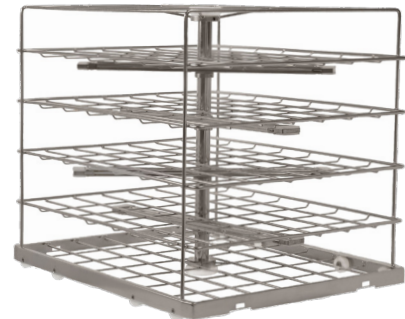


Fixed Rack

Six Levels

Ideal solution when large capacity processing is a request. Uses six washing levels.

- Washes up to 18 DIN instrument baskets per load.
- Maximum height of 2.3 inches (5.9 cm) between each level.

Item #P1500-FIX6

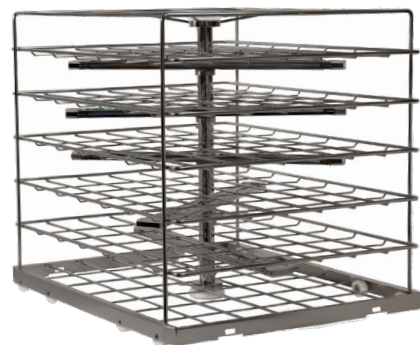


Anesthesia Rack

Tubes and Equipment

Designed to wash anesthesia tubes and equipment. Rack holds up to 12 anesthesia tubes, breathing bags, face masks, intubation instruments, laryngoscopes etc. 1 level suitable for DIN basket for face masks, intubation instruments, laryngoscopes and other equipment.

Item #P1500-AN

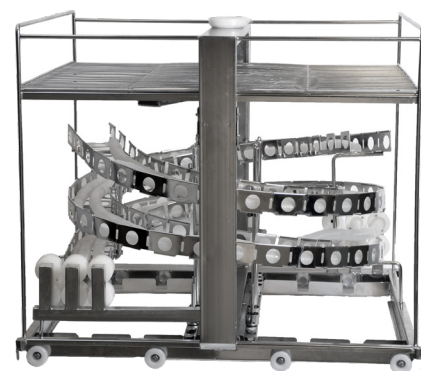

Wash Racks

Minimal Invasive Surgery Rack

instruments and rigid scopes

Rack for loading equipment guaranteeing efficient washing and disinfection of surgical instruments. Suitable for various OR-packages. The rack has room for e.g. equipment from two laparoscopic or biliary surgeries. Process up to 60 tubular instruments and 6 DIN baskets per load.

- 20 hoses
- 10 universal nozzles
- 10 nozzles with a long tip
- 20 nozzles for universal use

Item #P1500-MIS1



Minimal Invasive Surgery Rack

Instruments and rigid scopes

Wash racks are often used for Minimal Invasive Surgical Instruments, cleaning surgical instruments like scissors, forceps or retractors and special complex, hard-surfaced medical devices with narrow lumens such as suction cannulas, rigid endoscopes, drills, keys, wrenches, ratchets, insertion tools and hand held laparoscopic instruments.

- 8 scope spikes
- 8 holders for tubular instruments
- 1 holder for slim tubes
- 14 holders for suction and needles and similar
- 14 connections for suction hoses

Item #P1500-MIS2




Flexible Rack

Load carrier for vertical containers

The vertical container load carrier is designed to wash a number of containers vertically within the wash rack footprint. The design of the rack allows for movement of the retaining bars, making it easy to fit various sized containers.

- With a maximum height of 22 inches (55.9cm), the carrier can accommodate rigid containers of varying lengths

Item #P1500-CNTR

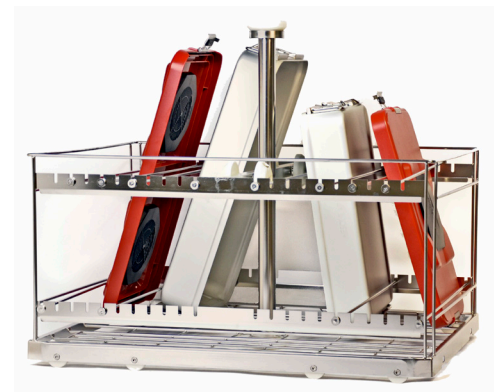

Rack Inserts

Insert Rack

Washing OR Shoes

Insert for efficiently washing up to 60 OR shoes per load.

- 60 pairs of hospital clogs or shoes
- Flexible and removable insert
- Used in conjunction with P1500-FLX2 two level wash racks



Item #P1500-FLX-SH

Shelf Level Rack

Washing Stainless Bowls

Supports the shelf divider arrangement for items of varying sizes.

- Washing capacity of 12 bowls of different sizes
- Maximum diameter of 14.2 inches (36 cm) per bowl
- To be combined with P1500-FLX2 and P1500-FLX4 wash racks



Item #P1500-B4L-BWL

Rack Inserts


Shelf Level Rack

Washing Plastic Bowls

Shelf level for Sterling P1500 wash racks for washing and disinfecting plastic bowls.

- Capable of holding up to 8 bowls
- Can hold bowls up to a max diameter of 13.8 inches (35 cm)
- To be combined with P1500-FLX2, P1500-FLX4, and P1500-FLX5 wash racks

Item #P1500-FLX-BWL

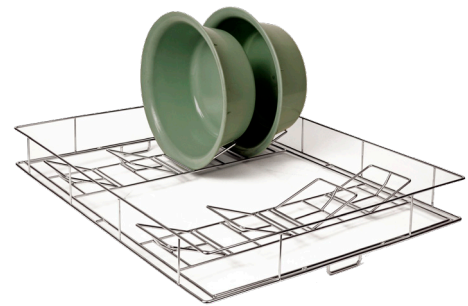


Shelf Level Rack

Washing Kidney Bowls

Shelf level for Sterling P1500 wash racks for washing and disinfecting kidney bowls.

- Capable of holding up to 18 kidney bowls
- 6 rails for various sized bowls
- To be combined with P1500-FLX2 and P1500-FLX4 wash racks

Item #P1500-FLX-KB

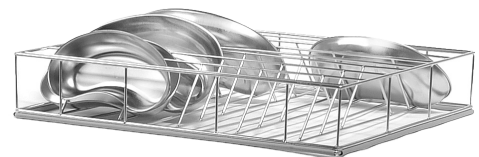



Circle of Infection Prevention Solutions



About Skytron: Skytron is a healthcare efficiency company specializing in capital equipment that drives performance in today's healthcare facilities. Our solutions enhance the utilization of people, facilities and capital because they are designed with the user in mind and have a low, long-term cost of ownership. We offer full-room solutions in Clinical, Infection Prevention and Clinical Business Intelligence products. **To learn more, visit our website at www.skytron.com**



5085 Corporate Exchange Blvd. SE
Grand Rapids, MI 49512
1.800.SKYTRON (759-8766)
email: info@skytron.com
www.skytron.com



REV-1 12-18
PN L2-010-138